




hp **hygie**[®]

ALU-CLASSIC
COMMODE CHAIR

**NEW
TECHNOLOGY**

“ **One** patient in twenty is impacted by a hospital-acquired infection ”

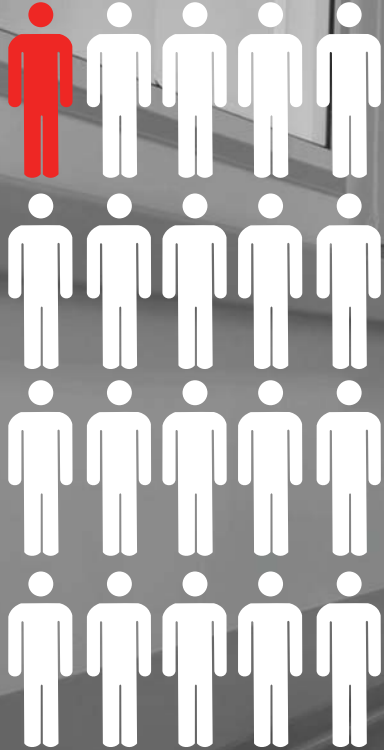

New technology in the prevention of infections

On the leading edge of the battle against nosocomial infections, Hygie is proud to introduce its brand new Alu-Classic® commode chair, made entirely of aluminium and free of porous surfaces. This Canadian-built product is specially designed for hospitals and long-term care facilities.

"This chair responds to infection control needs in health care facilities, as it contains no porous material. This helps facilitate cleaning and disinfection."

France - Infection Control Nurse

"The Alu-Classic commode chair is a must for hospitals and long-term care facilities. It's designed to facilitate access on both sides while providing patients with greater comfort and security."

Isabelle - Occupational therapist


"We are convinced that this new generation of commode chairs will significantly reduce patient falls. The rubberized wheels provide greater adherence and the four-wheel braking system ensures the chair remains completely motionless."

Claude - Nursing assistant



Technical specifications

①	Total width	57 cm	22.5 po
②	Total height	89 to 102 cm	35 in to 40 in
③	Total depth	68 cm	27 in
④ a	Floor to seat height (front)	38 to 43 cm	15 in to 17 in
④ b	Floor to seat height with extender (front)	46 to 51 cm	18 in to 20 in
⑤	Floor to seat height (back)	37 to 49 cm	14.5 in to 19,5 in
⑥	Seatback	42 cm	16.5 in
⑦	Seatback angle	5°	
⑧	Seat width	49 cm	19 in
⑨	Seat depth	44 cm	17.5 in
⑩	Distance between handles	15 cm	5.75 in
⑪	Clearance height	34 to 47 cm	13.5 in to 18.5 in
⑫	Arm rest height	20 cm	8 in
⑬	Wheel diameter	15 cm	6 in
⑭	Grip	13 cm	5 in
⑮	Side rail	51 cm	20 in
⑯	3 height adjustments	2.5 to 7.6 cm	1 in to 3 in
⑰	Floor to side rail height	20 cm	8 in
	Chair weight	9 kg	20,5 lb
	Maximum load	181 kg	400 lb

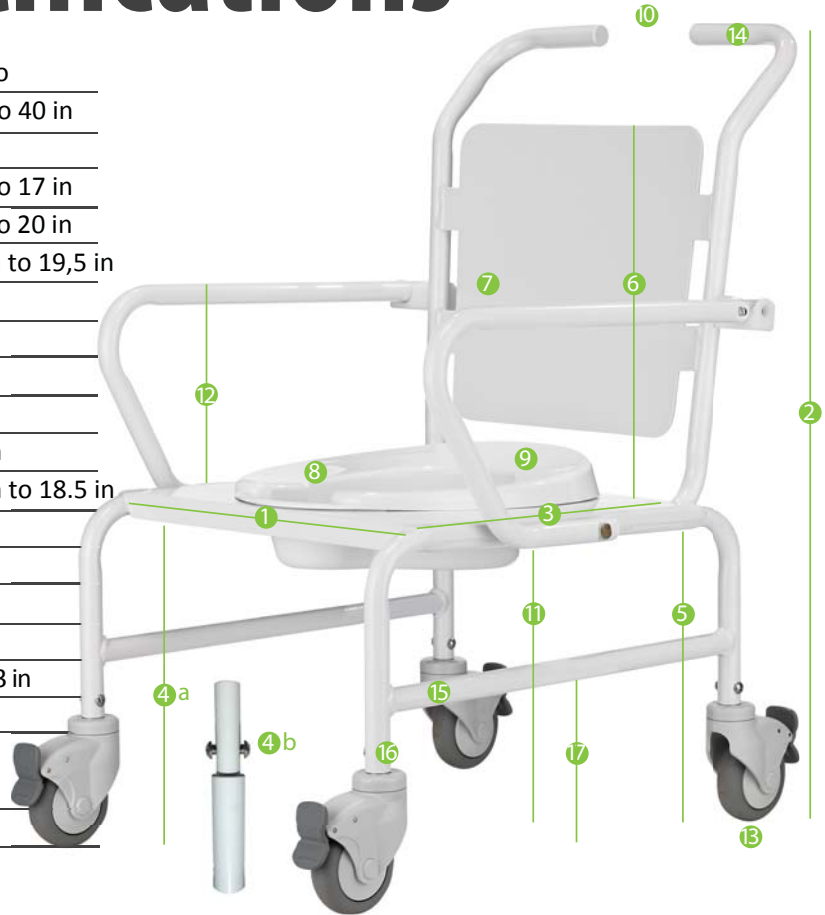

From 1990 to today

Over 80 options
to meet your
every need

HYGIE ALU-CLASSIC® CHAIR

- ✓ Aluminium construction
- ✓ Ultra-light, sturdy and stable
- ✓ Adjustable height
- ✓ Easy to clean
- ✓ Pivoting armrests
- ✓ Free of porous surfaces
- ✓ Rotates freely on itself
- ✓ Rust-free
- ✓ Floor to seat height of 15"



HYGIE CLASSIC® CHAIR

- ✓ Galvanized steel construction
- ✓ Sturdy and stable
- ✓ Adjustable height
- ✓ Easy to clean
- ✓ Pivoting armrests
- ✓ Free of porous surfaces
- ✓ Large size available
24" wide



WARRANTY: Chair frames are guaranteed for a period of two (2) years (Classic Chair) or five (5) years (Alu-Classic Chair) respectively from the date of purchase. Moving pieces (armrests, wheels, etc.) are guaranteed for a period of one (1) year from the date of purchase. Accessories (boots (armrest receptacles), basins, seats, etc.) are guaranteed for a period of three (3) months from the date of purchase. In the case of parts or pieces with a manufacturer's defect, Hygie® will supply replacement parts of original, similar or equivalent quality.

Wheelbase

The wheelbase of the Hygie Alu-Classic commode chair's frame makes it highly stable and secure for the patient, reducing the risk of overturning.




Rubber wheels

The 4-inch rubberized wheels offer 360° pivoting action and are equipped with individual manual brakes to fully immobilize the chair. The rubber also improves floor grip without leaving scuff marks.



Adjustable height

The extender allows for optimal, patient-tailored adjustments (see Technical Specs 4b). Three sizes available. This adjustment provides additional comfort for patients and a more comfortable posture for care staff.

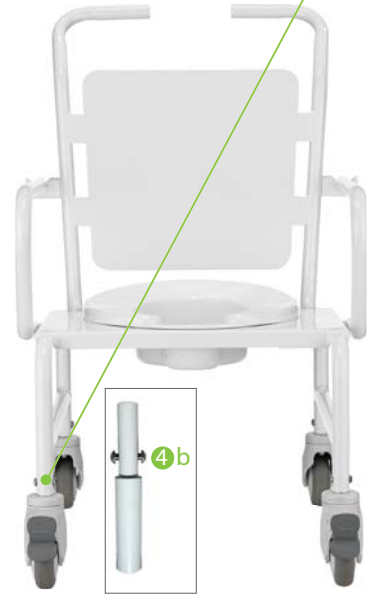

Pivoting armrests

The sturdy, secure armrests can be lowered by pressing a button. They also pivot along the seat axis, which frees that area and facilitates patient transfers.



Easy to clean

This chair was specially designed to minimize gaps and hard-to-reach spaces, ensuring quick and easy disinfection.



Construction

The lightness of the chair's frame is due to its aluminium construction. The slight tilt in the seatback provides better comfort to the patient.



Single-patient seats

Use of Hygie Hygienic Cover® recommended



AC-SIEG-HY21-000
Bedpan Seat with holes. Use with Hygie Hygienic Covers®

Cleaning bedpan seats and pails

- We recommend the use of detergents or enzyme products.
- Seats, basins and bedpans may be cleaned with a decontamination washer (maximum temperature of 90.5° C or 195° F).



1 CH-SIRO-OUSC-001
Round seat, open front without cover

2 CH-CUVE-HYAL-000
Pail with holes



AC-SIEG-HY21-SFO
Bottomless Bedpan Seat



AC-SIEG-CUVE-000
Bedpan Seat with bottom



Inside 1 Hygienic Cover
 Super-absorbent Pad | Absorbe and gelifie | 600ml
 Tampon super absorbant | Absorbe et gelifie | 30 sec.




 Made from recycled materials
 Fait de matières recyclées

Removable accessory plate with Hygienic Cover dispenser and a bracket to hold the disinfecting wipe receptacle. Install on the commode chair seatback.



Product code
CH-PLAQ-DIST-ALU

Hygienic Cover®
SA-HYGI-CLBI-000

The Hygie Hygienic Covers® for Commode Chairs are the perfect solution to effectively manage and contain body wastes at the point of care. They encourage healthcare providers to adopt optimal work practices, which reduces the spread of pathogens that cause epidemics, cross-contamination and nosocomial (hospital-acquired infections).

The super-absorbent pad in each envelope is unique – it can absorb more than 600 ml of body waste in seconds and transform it into a gel. Its use helps decrease the risk of splashes, spills and the spread of harmful germs, thus preventing nosocomial infections and cross-contamination.

Ease of maintenance and cleaning

- Designed to comply with new disinfection standards.
- To clean, use products with high detergent and disinfection properties such as a 5,000 ppm sodium hypochlorite solution following every use.
- Remove the bedpan prior to cleaning the frame.
- Begin the cleaning/disinfecting process, moving from least to most-soiled areas.
- We do not recommend the use of wheelchair washing machines. The mechanical action in manual cleaning is designed to enhance detergent efficiency by eliminating a high proportion of organic waste.
- Do not use or clean in the shower.
- Do not use a scouring pad.



Preventable statistics

10 % of patients admitted into a short-term care facility within the healthcare system will contract an infection as a result of received care or treatments.⁽¹⁾

30% to 50% of hospital-acquired infections are preventable.⁽²⁾

In Canada, **8,000 to 12,000** people die from nosocomial infections every year.⁽³⁾

Nosocomial infections cost the Quebec healthcare system approximately **\$180 million** annually.⁽⁴⁾

Nosocomial infections are the **fourth** leading cause of death in Canada.⁽⁵⁾

(1) <http://publications.msss.gouv.qc.ca/>

(2) Canadian Committee on Antibiotic Resistance, 2008.

(3) Zoutman, DE, Ford DB, Bryce E et al. The state of infection surveillance and control in Canadian Acute Care Hospitals; *Am J Infect Control*, 2003; 31:266-73.

(4) <http://publications.msss.gouv.qc.ca/>

(5) Statistics Canada. Ranking and number of deaths for the 10 leading causes, Canada, 2000 and 2009.






 1.866.588.2221

 www.hygie.com

 info@hygie.com

 9140 Leduc Blvd. Suite 230
Brossard (QC) Canada
J4Y 0L2



©2015, Hygie Canada inc. All rights reserved
Printed in Canada